




Scientific Equipment

Tel. 051 501153 - Fax 051 6336182

[www.favs.it](http://www.favs.it) - [info@favs.it](mailto:info@favs.it)

**IKA**

designed for scientists

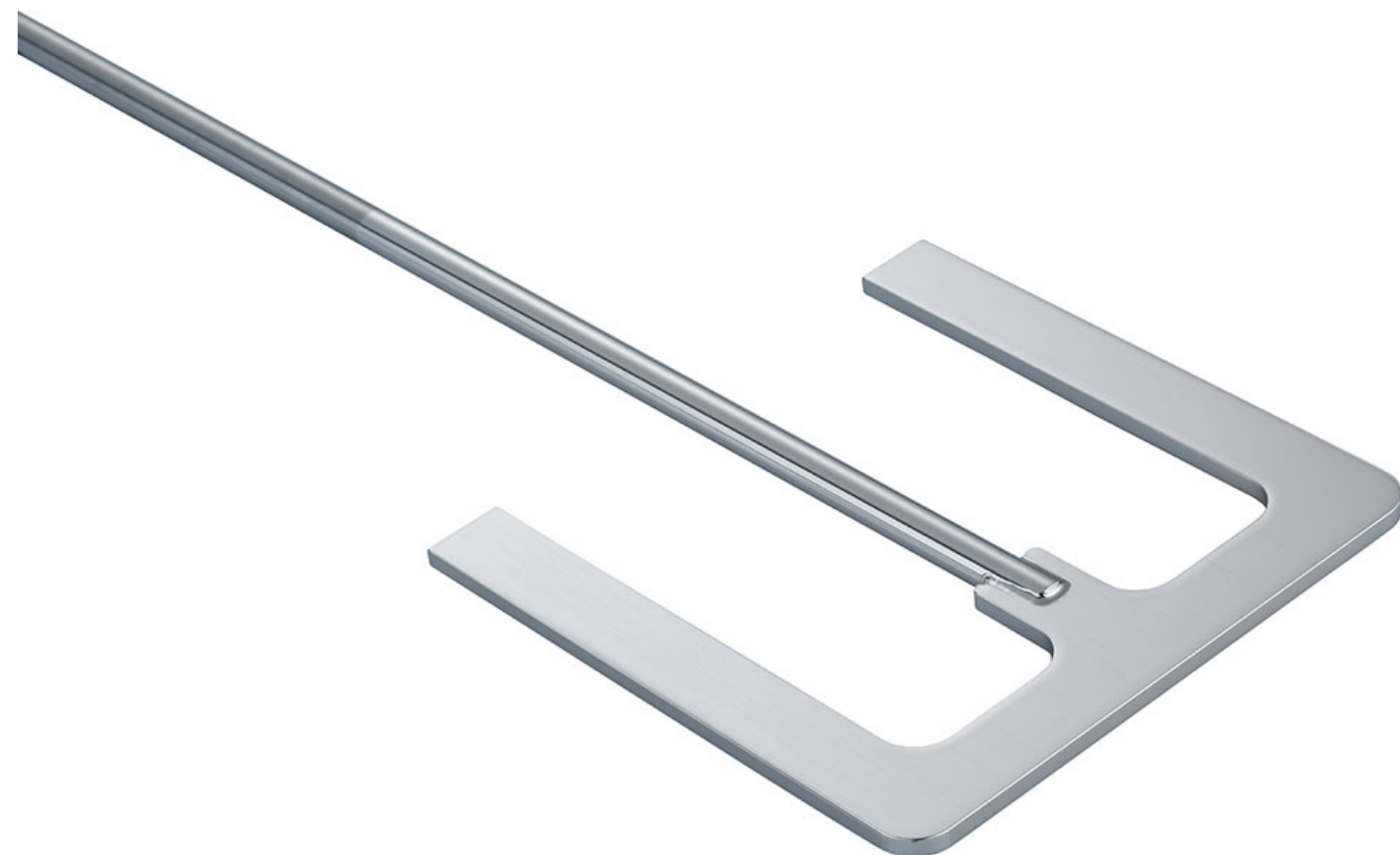

## R 1333 Anchor stirrer

/// Data Sheet

Tangential flow, high shearing rate at edges, minimum deposits on the vessel wall. Used at low speeds. Polymer reactions, even distribution of high mineral contents in liquids. The ideal stirrer for medium to highly viscous fluids. Stirrer

Ø: 150 mm

Shaft Ø: 10 mm

Shaft length: 550 mm

[www.ika.com](http://www.ika.com)

Subject to technical changes



IKAworlwide



IKAworlwide /// #lookattheblue



@IKAworlwide



designed for scientists

### Technical Data

Stirrer diameter [mm]	150
Shaft diameter [mm]	10
Shaft length [mm]	550
Speed max. [rpm]	800
Material in contact with medium	AISI 316L / 1.4571
Weight [kg]	0.638

